

Logging SOP Assessment 1

This quiz is open note/open wiki, but it must be taken individually. Please do not share information about the quiz contents with your coworkers.

The SPS Wiki can be found here: [Lab Med Wiki](#)

Question 1.

Have you read all 8 Reference Materials and all 10 General Logging SOPs since they were updated on February 1st, 2013?

- Yes
 No

Question 2. [Question Pulled From: SP.PR.S.AZ Shared Specimens and Shared Results](#) [Now Also Listed Under: SP.PR.S.AZ Test Cancellation](#)

You accidentally routed a BMP to a lime tube instead of a gold one, so you log it to a new accession number. The original collection time was 13:02. What collection time do you enter in the new accession number?

- 13:02
 13:03
 UNKNOWN
 13:05

Question 3. [Question Pulled From: SP.PR.S.AZ Routing Containers and Resulting on Entry](#)

One urine container is received for a Drug Screen and a Urine Pregnancy. After receiving the order in ORM, the Misys Routing Box defaults to two random urine (UR) containers. What do you do?

- Press 'Route' and place all labels on the urine container.
 Press 'Route', print an extra set of labels, and place all labels on the urine container.
 Remove one UR and route testing onto only one container. Aliquot the container and place the additional labels onto the specimen for the Line Loader to deliver.
 Remove one UR and route testing onto only one container. The Spin technician will aliquot the container.

Question 4. [Question Pulled From: SP.PR.S.AZ Creating and Selecting Episodes / Encounters](#)

You receive a requisition with patient information but no encounter number or location. What do you do...


for an inpatient?

Call Admitting.

for an outpatient?

Give to a Lead. The lead will find the patient location in Epic.

Question 5. [Question Pulled From: SP.PR.S.AZ Routing Containers and Resulting on Entry](#)

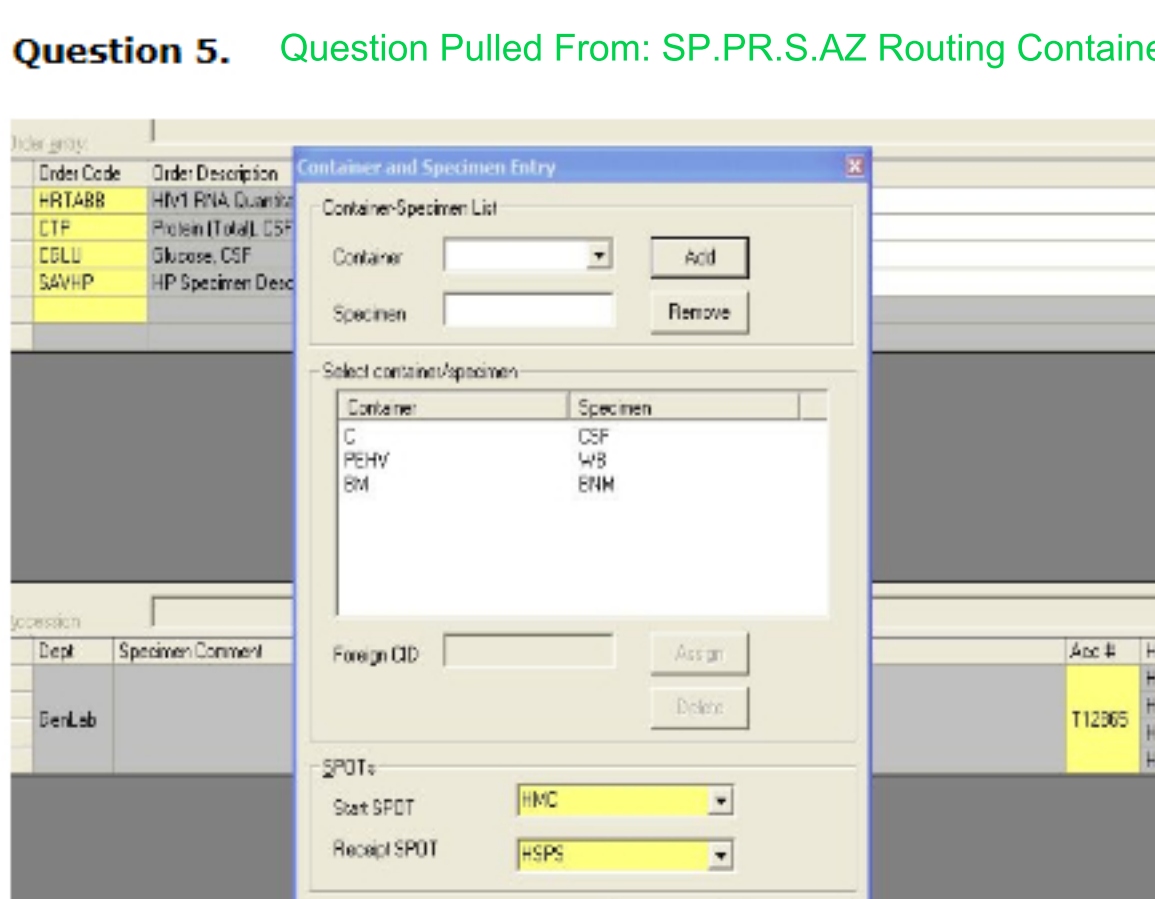

You receive a syringe of CSF for the following tests: HRTABB, CGLU, and SAVHP. After you save the order, the Specimen and Container Entry box appears with the following containers defaulted: C, PEHV, BM. What containers should you route?

- C,PEHV,SP
 C,C,SP
 C
 C,PEHV,BM
 C,C,C
 C,SP
 C,C

Question 6. [Question Pulled From: SP.PR.S.AZ Shared Specimens and Shared Results](#)

An outpatient arrives at the lab for a blood draw with three orders:

1. CBC ordered by Dr. Dumbledore from the Hogwarts Clinic,
2. CBD ordered by Dr. Pevensie from the Narnia Clinic,
3. CBC ordered on a green RTS from the Cat Flu Study.

After creating the order, what additional steps must be taken? (Check all that apply)

- Mark the CBC on order #1 to be credited
 Mark the CBC on order #3 to be credited
 Fill out a QA on the Specimen Tracker
 Cancel and credit the CBC on order #3 as a duplicate
 Fill out a Shared Specimen/Results form
 Mark the CBC on order #2 to be credited
 Cancel and credit the CBC on order #1 as a duplicate

Question 7. [Question Pulled From: SP.PC.S.REF Cancel Codes](#)

You receive a blue top for a CBC order. You call RN Arnie at 1234 and are told the order will be recollected.

What cancel code should you use?

- YIMP
 YNSR
 YPROV
 YCOLL

Question 8. [Question Pulled From: SP.PR.S.AZ Test Cancellation](#)

What free-text message should be entered after the cancel code?

- ;ARNIE
 ;RN ARNIE WAS NOTIFIED
 ;RN Arnie notified at 1234
 ;RN Arnie will recollect

Question 9. [Question Pulled From: SP.PC.S.REF Cancel Codes](#)

Select the correct Cancel Code per the Cancel Code SOP.

Order changed by provider

after receipt in SPS

(includes cancelation due to

reschedule).

YORDC

Specimen integrity

compromised in SPS.

YLACC

Specimen was mislabeled in

the lab and results were

reported.

YDRWC

Order canceled by provider

before receipt in SPS

(includes cancelation due to

reschedule).

YPROV

Question 10. [Question Pulled From: SP.PC.S.REF Master Forms](#)

Go to the Master Forms SOP. Open the "Patient Instructions for Timed Urine". What is the first tip listed under "Tips for successful collection"?

Note: Instructions are the same at both locations.

- Collect all of the urine your body produces during the 24 hour period.
 After 24-hour collection, please bring full collection bottle back to the laboratory for testing
 Avoid alcoholic beverages and vitamins for 24 hours before and during the collection of your urine.
 During collection, 24-hour urine collection bottle should be refrigerated
 Do not discontinue taking medications during your collection unless directed by your healthcare provider or
 Begin and end your urine collection at the same time, 24 hours apart.

Question 11. [Question Pulled From: SP.PC.S.REF Useful Logging References](#)

According to the Logging References SOP, which one of these handwritten requests requires clarification?

- LFTs
 BAL
 Cardiac Enzymes
 Nutrition Labs
 CMP

Question 12. [Question Pulled From: Section 1. Reference Materials \(SP.PC.S.REF\)](#)

Where can you find each of the following:

A Minutes-to-Decimals

Conversion Table

SP.PC.S.REF Decanting References

A list of test requests that

require clarification

SP.PC.S.REF Useful Logging References

Jenn Reichert's phone

number

SP.PC.S.REF Useful Phone Numbers

The routing code for a dark

green tube

SP.PC.S.REF Tube Types and Routing Codes

Question 13. [Question Pulled From: SP.PR.S.AZ Routing Containers and Resulting on Entry](#)

A COMP, HFPA, MG, PPP, and CBC are requested on the same requisition. The PPP and HFPA appear on accession number **W34567**, the COMP and CBC are on accession number **H87654**, and the MG is on accession number **F12345**. You received one Blue and one Lime tube. Which is the most correct way to route the containers?

Requisition: COMP, HFPA, MG, PPP, CBC

Tubes: 1 BL tube, 1 PS tube

Accession Numbers:

W34567 - PPP, HFPA

H87654 - COMP, CBC

F12345 - MG

- Remove PS and add NT to both W34567 and H87654 before canceling all tests as YDUP. Add a COMPHF to F12345.
 Receive and cancel all three accession numbers as they are. Cancel all tests as YCPFIX. Then, create a new order.
 Remove the PS and add NT on both W34567 and F12345. Cancel the MG and HFPA as YCPFIX. Then add BILD and MG to H87654. Remove the L tube from H87654 and route to NT. Cancel the CBC as YDATA.
 Remove the PS on W34567 and route only the BL tube. Also Remove the PS on F12345, but route as NT. Cancel the MG and HFPA as YCPFIX. Then, add BILD and MG to H87654. Remove the L tube from H87654 and route to only one Lime tube. Cancel the CBC as YIMP.