
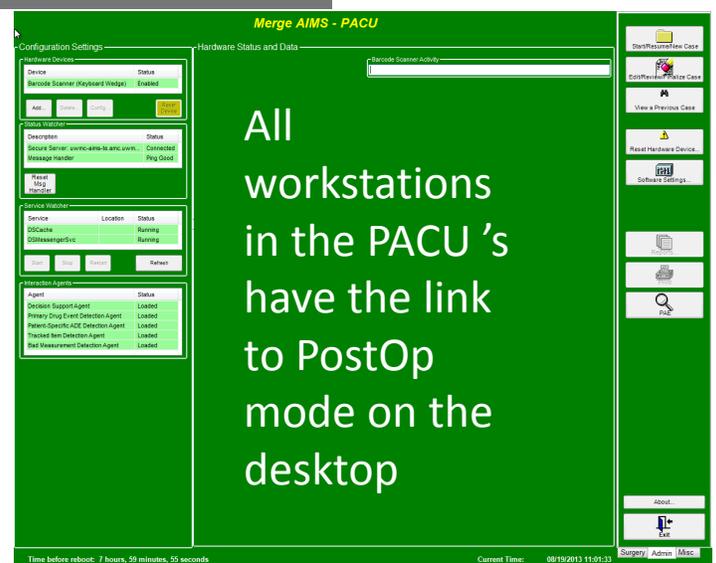
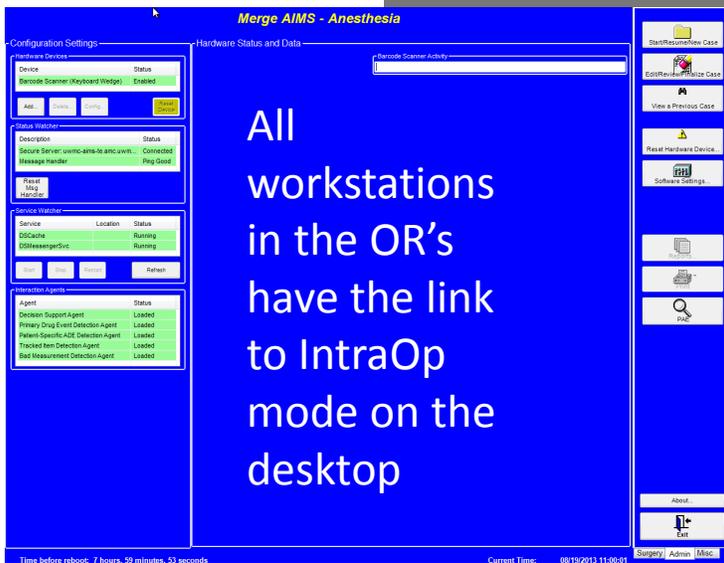
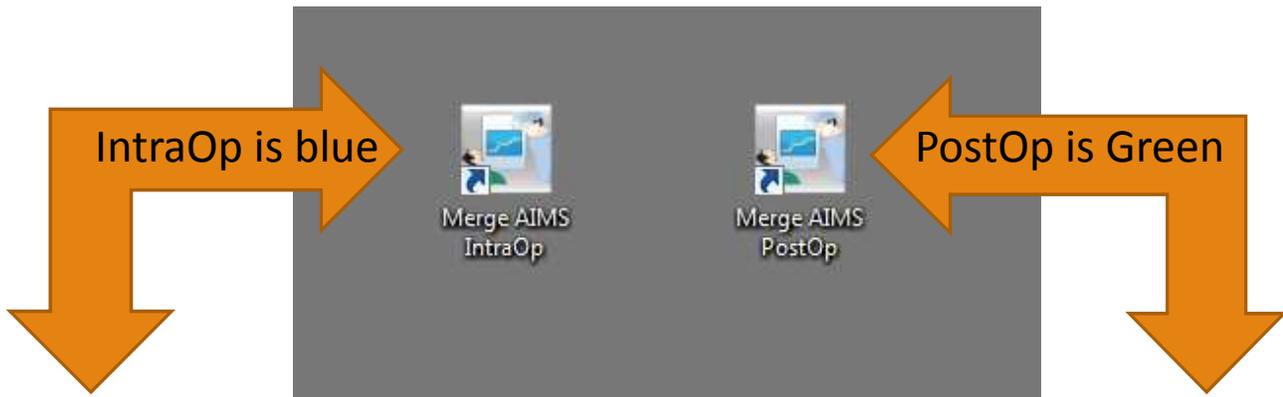


Merge Tutorial

The main desktop layout of all Merge AIMS computers have a grey background to help differentiate them from standard desktops.

There are two types of workstation modes one is OR and the other is PACU. Both are still grey in color but have different icon links. (**IntraOp/PostOp** mode)

IntraOp mode will be a blue background and PostOp mode will be a green.

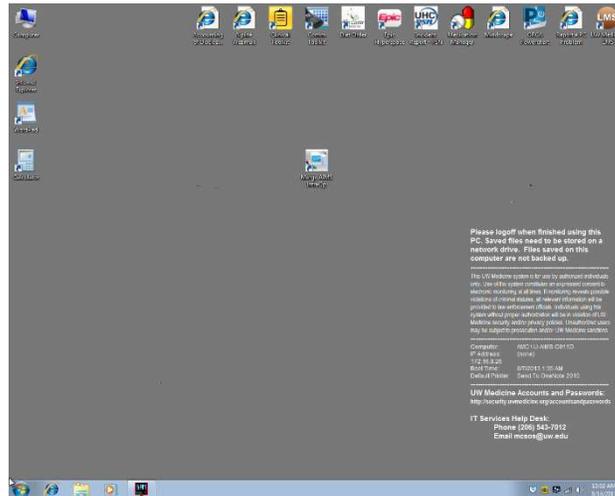

Starting a case

Its very important to make sure all system checks are working

Configuration Settings

Hardware Devices

Device	Status
Monitor (DataCaptor - localhost:7272)	Enabled

Buttons: Add..., Delete..., Config..., Reset Device

Status Watcher

Description	Status
Secure Server: uwmc-aims-ap.amc.uwm...	Connected
Monitor (DataCaptor - localhost:7272)	Connected
Message Handler	Ping Good

Buttons: Reset Msg Handler

Service Watcher

Service	Location	Status
DSCache		Running
DSMessengerSvc		Running

Buttons: Start, Stop, Restart, Refresh

Interaction Agents

Agent	Status
Decision Support Agent	Loaded
Primary Drug Event Detection Agent	Loaded
Patient-Specific ADE Detection Agent	Loaded
Tracked Item Detection Agent	Loaded
Bad Measurement Detection Agent	Loaded

Merge AIMS - Anesthesia

Hardware Status and Data

Record Map	Measurement	Value	Time	Source
------------	-------------	-------	------	--------

Buttons: Refresh Data, Show Hidden Measurements Only

Time before reboot: 7 hours, 39 minutes, 54 seconds
Current Time: 08/07/2013 09:14:25

Please verify all Status indicators are green before starting a case to insure proper communication between AIMS applications servers.




Any other color RED or YELLOW Please notify Aims Support before starting a case.

PH: 206.598.8695

PGR: 206.540.0937

All Items in the above image need to be green before starting a case. This is on the main merge startup page

Starting a case

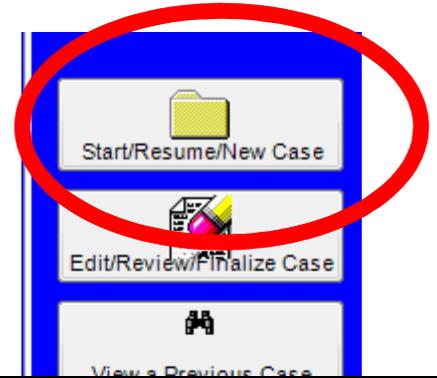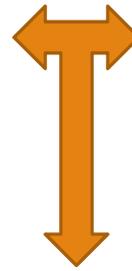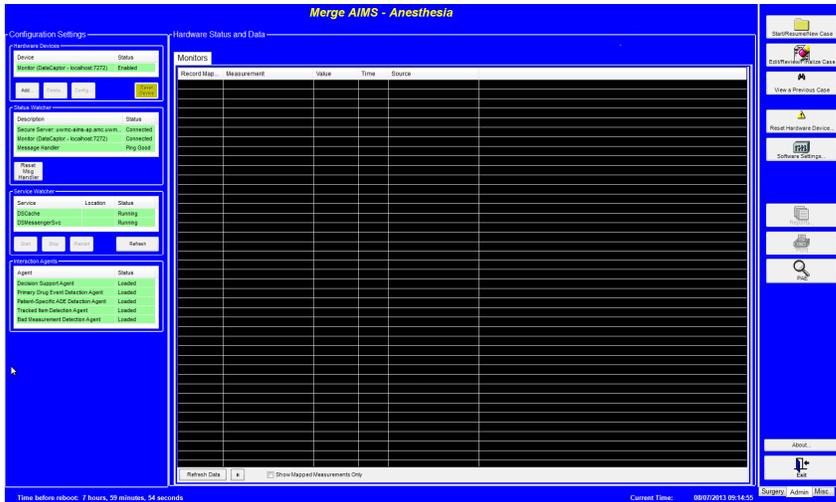

Login required for Authenticated Case Close

User Name:

Password:

If all indicators are green start a case by pressing the Start/Resume/NewCase button. The system will ask you to login (AMC-User name and Password) same as ORCA login.

Start, Resume or Create a Case

Patient	Procedure	Time	Acct. Num.	MRN	Loc.
ZZ, BUBBA	CABG & MVR	08:55 08/14/2013	1763274340	U3210037	Test OR
ZZ, CUTEBUNNY	CYSTOSCOPY W/...	09:00 08/15/2013	1763277571	U3253483	Test OR
ZZ, PENNY	CABG W/ IMA	08:30 08/14/2013	1763275504	U3016162	Test OR
ZZSTORM, SEATTLE	CYSTOSCOPY, U...	07:30 08/15/2013	1763277570	U4006766	Test OR
ZZTEST, MAGA	APPENDECTOMY ...	07:00 08/15/2013	U3000001	U3000001	Test OR


List Filtering:
 Show for Date Range Exclude Closed Cases
 Show for Staff: VERGARA, R
 Show for Patient
 Show for Location: Test OR

Date Range:
Starting Date: 08/12/2013
Ending Date: 08/19/2013

Last Name:
Search String (Partial OK)

Selected Procedure: ZZ, CUTEBUNNY (Acct #: 1763277571, DOB: 01/01/1951), CYSTOSCOPY W/ BLADDER BIOPSY 8/15/2013 9:00:00 AM. Case not Started.

If your patient is not on this list try clicking the update procedure button. If still not on the list select "Emergency /Unkown Patient. (See instruction on how to Identify TBD case.)

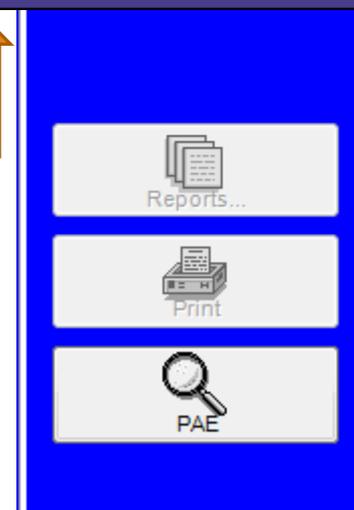

When logged in select the desired patient from the case list and press OK and confirm your selection before proceeding.

Confirm Case Selection

Are you sure this is the case you intended to select?

Patient: ZZ, BUBBA KIMBERLEE
Acct #: 1763274340, MRN: U3210037
DOB: 08/08/1880, SSN:

Procedure: CABG MVR (CABG MVR*)
Procedure Time: 08/14/2013 08:55